

DUOBLADE

ADVANCED ENERGY BLADE

FOR MINIMAL THERMAL INJURY,
LESS BLEEDING, REDUCED SURGICAL SMOKE



A D V A N C E D E N E R G Y B L A D E

DUOBLADE™ DEVICES

DUOBLADE™ AIMS TO REDUCE SURGICAL RISKS AND ENHANCE EFFICIENCY

Lower Temperature

With a self-limiting mechanism, the blade stays below 100°C, reducing the operating temperature by over 60% compared to traditional electrosurgery.

Better Hemostasis

Providing more efficient surgery by minimizing thermal damage while maximizing the hemostatic performance inherent to electrosurgery.

Incision Healing

Markedly minimizes the depth of thermal injury and the width of scars in comparison with traditional electrosurgery.

Less Surgical smoke

Remarkably reduced smoke, ensuring better visibility and protecting the health of medical staff.

Product model	Description
DB1SEP	Duoblade™ with smoke evacuation, rotating and telescoping shaft
DB1SP	Duoblade™ with rotating and telescoping shaft
DB1EP	Duoblade™ with smoke evacuation
DB1P	Duoblade™ essence

Blade selection	Description
Standard tip	Standard Blade Electrode
Distal tip	Only 3mm of distal tip exposure which may protect alternate site thermal injuries.



「DUO」 BLADE™ TECHNOLOGY

Dielectric **U**ltra-focused **O**scillatory Heating


Minimal thermal damage



DUO TECHNOLOGY FOR FASTER RECOVERY & SAFE SURGERIES

Duoblade™ utilizes Dielectric Ultra-focused Oscillation (DUO) technology, a combination of a uniquely shaped cutting-edge blade fully insulated with a heat-resistant oxide layer and non-stick coating.


This localizes a uniform electric field around the edge of the blade with a concentration at the distal tip, enabling exceptional low temperature soft tissue dissection blended with simultaneous optimized hemostasis around the blade perimeter.

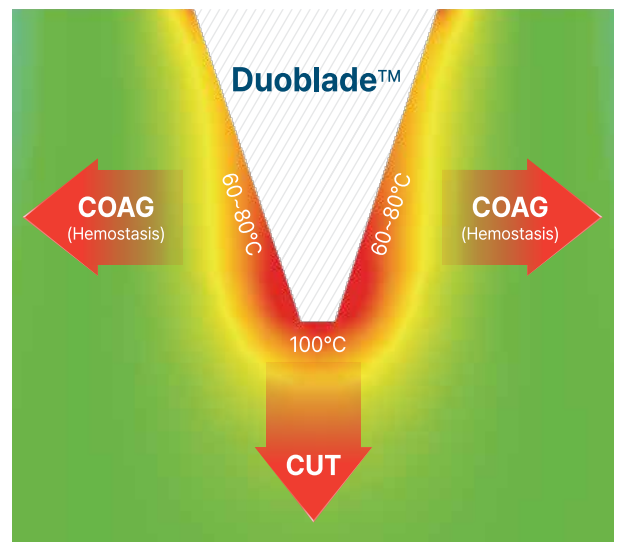

Electric field distribution and temperature effect of Duoblade™

THE LOWEST CUTTING OPERATING TEMPERATURE¹

The Duoblade™ features a uniquely shaped cutting edge that is doubly insulated with an electrochemically converted, stable, and heat-resistant oxide layer. This advanced technology ensures exceptionally low-temperature soft tissue dissection, automatically limiting the tissue temperature to below 100°C. These infrared thermal images show that the traditional electrosurgery device has a high operating temperature.

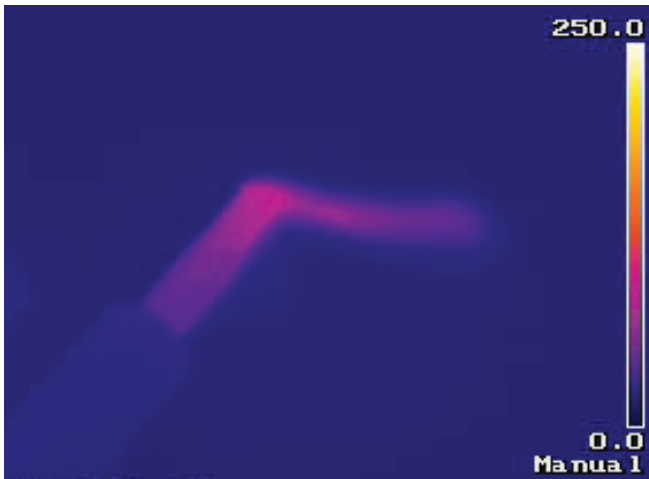


Duoblade™ operates up to 100°C




Traditional electrosurgery devices : 230 – over 280+°C

HISTOLOGY IMAGE PROFILE¹

Up to 70% reduction in thermal injury compared to conventional energy devices. These histologic profiles compare the thermal injury (marked in red) between the Duoblade™ and traditional electrosurgery devices under similar cutting settings.

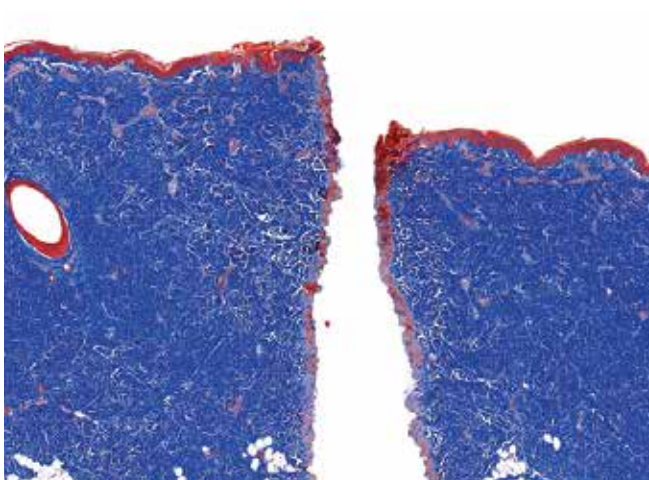

Duoblade™ incision using the CUT 50W setting, showing low thermal injury.

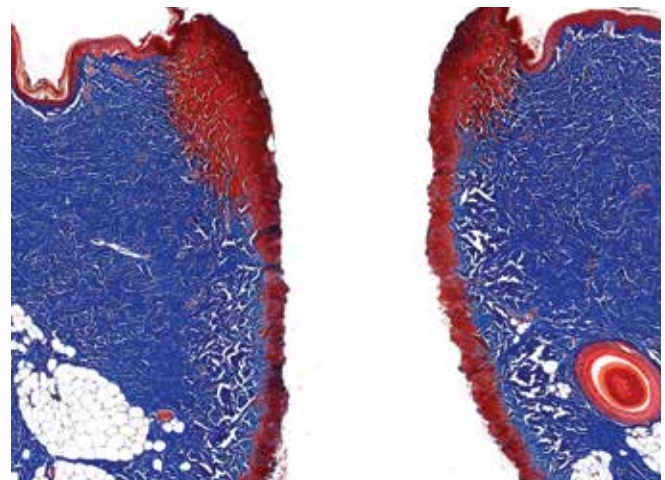


Incision using a traditional electrosurgical tool using the CUT 50W setting, showing significant thermal injury.

¹ TRE-00069 Rev. A. Preclinical operating temperature study. Data on File


SUPERIOR HEMOSTATIC PERFORMANCE²

Throughout the quantification of blood loss, outstanding hemostatic characteristics were shown at the CUT & COAG 35W settings. Duoblade™ demonstrated significantly smaller estimated blood loss compared to a scalpel and other competitive devices (P < 0.05).

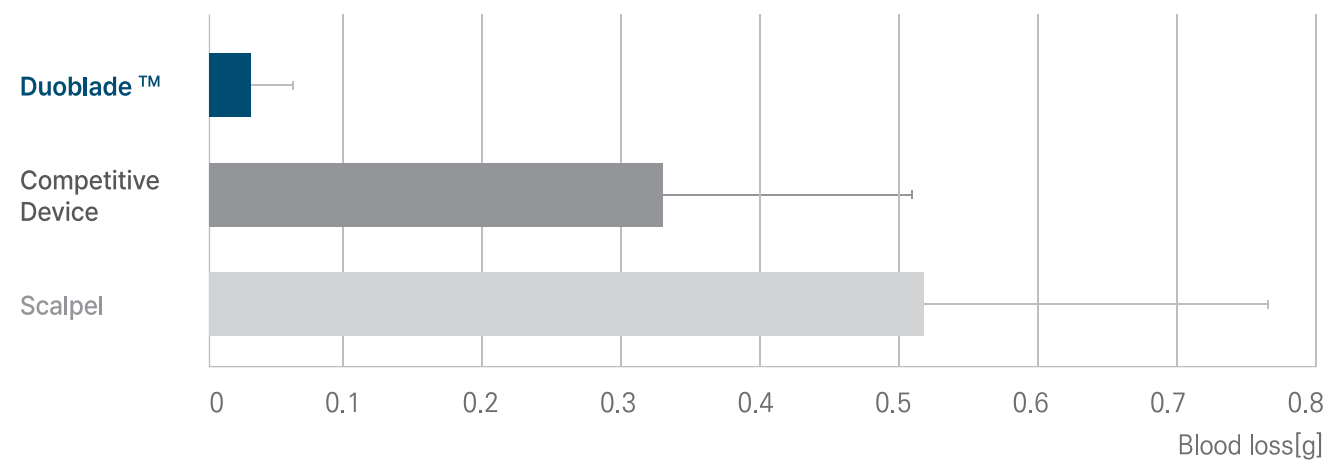

MINIMAL SURGICAL SMOKE³

Minimized surgical smoke due to the low operating temperature, reducing inhalation risks for OR staff through optimized-temperature operation. Air around the incision site (60 liters) was collected, and the number of particles was measured for surgical smoke sizes of 0.3 µm and 0.5 µm.

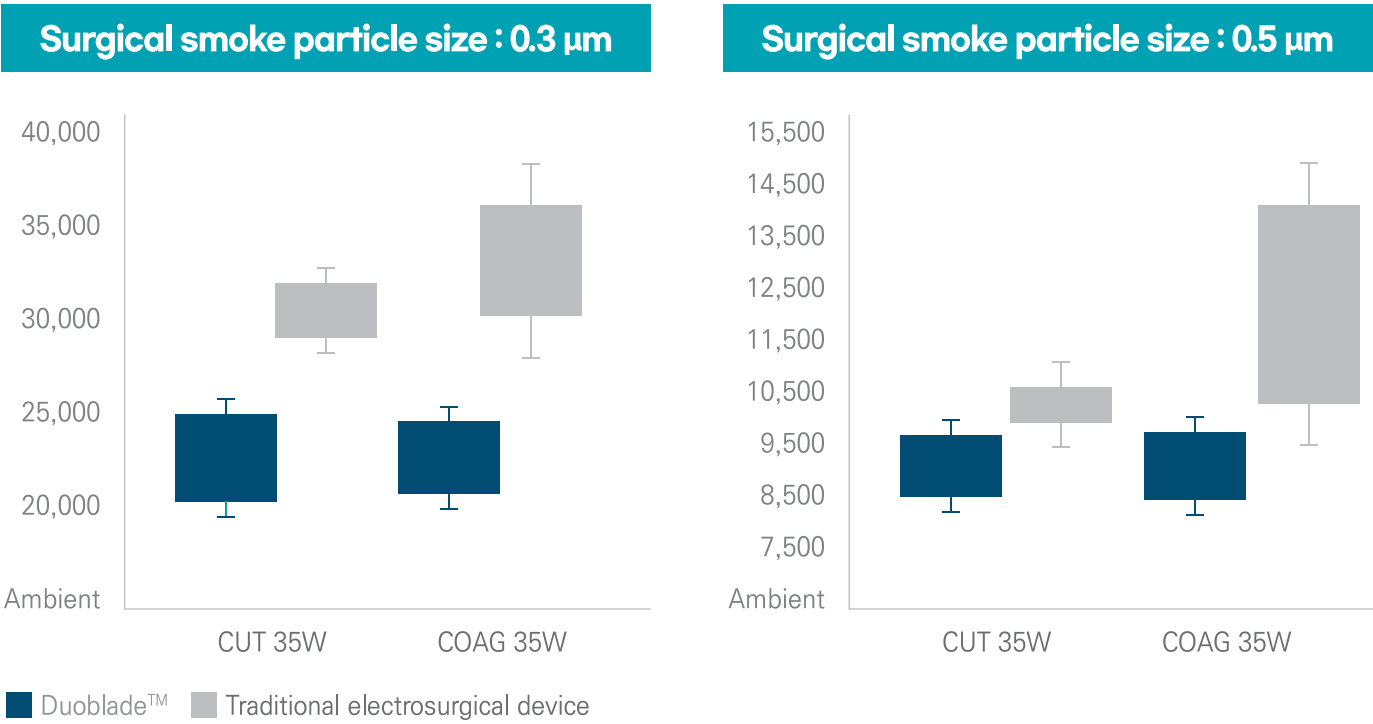

² TRE-00082 Rev. A Comparative Bleeding Control Study. Data on File
³ DOC-00044 Duoblade Technical Deck

© CRESEN 2025–2026. All rights reserved.

This document is intended specifically for healthcare professionals. Surgeons must rely on their own professional clinical judgment when deciding whether to use a particular product for treating a specific patient. We do not provide medical advice and recommend that surgeons receive proper training in the use of any product before employing it in surgery.

The information provided is intended to illustrate CRESEN's products. Surgeons should always refer to the package insert, product label, and instructions for use, including any cleaning and sterilization guidelines (if applicable), before utilizing any of CRESEN's products.

Product availability may vary by market due to differing regulatory and medical practices. Please contact your representative if you have questions regarding the availability of CRESEN's products in your region.

Reproduction, in whole or in part, is prohibited without the written consent of the copyright holder. Specifications and features are subject to change without notice. All changes will comply with regulations governing the manufacture of medical equipment.

CRESEN owns, uses, or has applied for trademarks or service marks, including Duoblade and CRESEN. The absence of a product, feature, service name, or logo from this list does not constitute a waiver of CRESEN's trademark or other intellectual property rights concerning that name or logo.



6200 Stoneridge Mall Rd. 3rd Fl.
Pleasanton, CA 94588 USA

www.cresen.co
hello@cresen.co
619-356-0583



MKT-00084 Rev. B